ALFA 6

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GRAZIE!

DISPOSAL OF THE ELECTRIC OR ELECTRONIC DEVICE AT THE END OF ITS LIFE

**INFORMATION FOR USERS ABOUT THE CORRECT DISPOSAL
OF THE PRODUCT ACCORDING TO THE EUROPEAN DIRECTIVE 2012/19/EU
AND TO THE LEGISLATIVE DECREE N. 49/2014.**

DEAR CUSTOMER,

IT IS IMPORTANT TO REMIND YOU THAT, AS A GOOD CITIZEN, YOU CANNOT DISPOSE OF THIS PRODUCT AS MUNICIPAL WASTE; YOU HAVE TO DO A SEPARATE COLLECTION. THE SYMBOL LABELLED ON THE APPLIANCE INDICATES THAT THE RUBBISH IS SUBJECT TO "SEPARATE COLLECTION". ELECTRICAL AND ELECTRONIC APPLIANCE CANNOT BE CONSIDERED AS DOMESTIC WASTE. THEREFORE, AT THE END OF ITS USEFUL LIFE, THIS PRODUCT SHOULD NOT BE DISPOSED WITH MUNICIPAL WASTES. THE USER MUST THEREFORE ASSIGN OR (HAVE COLLECTED) THE RUBBISH TO A TREATMENT FACILITY UNDER LOCAL ADMINISTRATION, OR HAND IT OVER TO THE RESELLER IN EXCHANGE FOR AN EQUIVALENT NEW PRODUCT.

THE SEPARATE COLLECTION OF THE RUBBISH AND THE SUCCESSIVE TREATMENT, RECYCLING AND DISPOSAL OPERATIONS PROMOTE PRODUCTION OF APPLIANCES MADE WITH RECYCLED MATERIALS AND REDUCE NEGATIVE EFFECTS ON HEALTH AND THE ENVIRONMENT CAUSED BY IMPROPER TREATMENT OF RUBBISH.

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EN USER INSTRUCTIONS

FOR A CORRECT AND SAFE USE OF THE APPLIANCE

K O O P E R[®] X

ONLY FOR A DOMESTIC USE

ATTENTION! READ THE INSTRUCTIONS GIVEN IN THIS HANDBOOK VERY CAREFULLY.

IT SHOULD BE PUT ASIDE FOR FURTHER CONSULTATION DURING THE ENTIRE LIFE-CYCLE OF THE PRODUCT SINCE IT GIVES IMPORTANT INDICATIONS ON THE PRODUCT'S FUNCTIONS AND SAFETY DURING INSTALLATION, USE AND MAINTENANCE.

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GENERAL SAFETY INDICATIONS

~~DANGER FOR CHILDREN~~ IS NOT A TOY. AFTER USE, STORE IT OUT OF REACH OF CHILDREN. KEEP PACKING MATERIAL (PLASTIC BAGS, CARTON, ETC.) OUT OF REACH OF CHILDREN.

- THIS APPLIANCE IS INTENDED TO BE USED BY CHILDREN 8 YEARS OF AGE OR OLDER AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, SUPERVISED OR INSTRUCTED BY ADULTS RESPONSIBLE FOR THEIR SAFETY. THIS APPLIANCE IS NOT INTENDED FOR USE BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, UNLESS THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY AND HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY.
- CHILDREN SHOULD BE SUPERVISED TO MAKE SURE THEY DO NOT PLAY WITH THE APPLIANCE.

DANGER OF INJURY!

- PAY UTMOST ATTENTION IN CASE THE APPLIANCE IS USED IN THE PRESENCE OF CHILDREN.
- NEVER LEAVE THE APPLIANCE UNATTENDED WHILE IN USE.
- TO AVOID INJURIES, KEEP CLOTHES AND ANY OTHER OBJECTS AWAY.
- SWITCH OFF AND UNPLUG THE APPLIANCE BEFORE CHANGING THE ACCESSORIES! DO NOT PLUG THE APPLIANCE TO THE POWER SOCKET IF THE NECESSARY ACCESSORIES ARE NOT INSTALLED.
- NEVER TOUCH THE MOVING PARTS OF THE APPLIANCE AND ALWAYS WAIT FOR THE APPLIANCE TO STOP COMPLETELY.
- UNPLUG THE APPLIANCE FROM THE POWER SOCKET WHEN THE APPLIANCE IS NO LONGER IN USE.
- DO NOT CLEAN WHILE THE APPLIANCE IS WORKING.
- THE CLEANING MUST NOT BE CARRIED OUT BY UNSUPERVISED CHILDREN.

DANGER OF ELECTRICAL SHOCK!

- THIS APPLIANCE IS DESIGNED TO OPERATE ONLY WITH ALTERNATING CURRENT. BEFORE STARTING IT, MAKE SURE THAT THE VOLTAGE CORRESPONDS TO THE ONE INDICATED ON THE RATING LABEL.
- CONNECT THE APPLIANCE ONLY TO THE POWER OUTLETS WITH A MINIMUM CAPACITY OF 10 AMPERES.
- DO NOT PULL THE POWER CORD OR THE APPLIANCE ITSELF TO REMOVE THE PLUG FROM THE WALL SOCKET.

- FULLY UNWIND THE POWER CORD BEFORE OPERATING THE APPLIANCE.
 - NEVER LET THE POWER CORD RUN THROUGH SHARP EDGES; AVOID CRUSHING IT.
 - DO NOT LEAVE THE APPLIANCE EXPOSED TO ATMOSPHERIC AGENTS (RAIN, FROST, SUN, ETC.).
 - NEVER USE THE APPLIANCE:
 - WITH WET OR MOIST HANDS;
 - IF IT HAS FALLEN AND HAS OBVIOUS RUPTURES OR MALFUNCTIONS;
 - IF THE PLUG OR POWER CORD DEFECTIVE;
 - IN HUMID OR GAS SATURATED ENVIRONMENTS.
 - DO NOT COVER THE APPLIANCE WITH TOWELS, CLOTHES, ETC.
 - PLACE THE APPLIANCE ONLY ON HORIZONTAL AND STABLE SURFACES, AWAY FROM SHOWERS, BATHTUBS, SINKS, ETC.
 - IF THE APPLIANCE GETS WET OF GETS IN CONTACT WITH WATER, DO NOT USE IT AND MAKE IT CHECKED BY A SPECIALISED LABORATORY.
 - ALWAYS REMOVE THE MAINS PLUG FROM THE POWER SUPPLY IMMEDIATELY AFTER USE AND BEFORE PERFORMING ANY CLEANING OR MAINTENANCE WORK.
 - IF POWER CORD IS DAMAGED, IT HAS TO BE REPLACED WITH ANOTHER CABLE OR A SPECIAL ASSEMBLY BY THE MANUFACTURER OR BY A CUSTOMER SERVICE OR BY SOMEONE WITH A SIMILAR QUALIFICATION, IN ORDER TO AVOID ANY RISK.
 - REPAIRS SHOULD BE CARRIED OUT ONLY BY A QUALIFIED TECHNICIAN.
- PLEASE NOTE THAT REPAIRING THE APPLIANCE BY

YOURSELF CAN BE EXTREMELY DANGEROUS.

- AN APPLIANCE THAT IS NOT WORKING ANYMORE, MUST BE RENDERED UNUSABLE BY CUTTING THE POWER CORD BEFORE DISPOSING OF IT.
- CLEAN THE APPLIANCE USING A DAMP CLOTH. DO NOT USE SOLVENTS OR ABRASIVE DETERGENTS.
- THE APPLIANCE IS DESIGNED FOR A DOMESTIC USE.

ATTENTION!

THIS APPLIANCE IS NOT A SEX TOY: DO NOT USE IT ON INTIMATE PARTS AND / OR GENITALS.

THIS APPLIANCE IS NOT A MEDICAL DEVICE AND THUS IT MUST NOT BE USED NEITHER FOR THERAPEUTIC OR SIMILAR PURPOSES NOR A PROFESSIONAL USE.

PLEASE READ THIS INSTRUCTION MANUAL CAREFULLY BEFORE INSTALLING AND USING THE APPLIANCE. KEEP IT FOR FUTURE REFERENCE.

BATTERY WARNINGS

- DO NOT IMMERSE ANY PART INTO WATER.
- AND DO NOT EXPOSE TO HEAT OR FIRE.

DANGER OF EXPLOSION!

- DO NOT PIERCE OR TAMPER WITH THE BATTERIES;
 - DO NOT THROW THE BATTERIES INTO THE FIRE;
 - IF THE LEAKING FLUID FROM BATTERIES GETS IN CONTACT WITH CLOTHES OR SKIN, WASH IMMEDIATELY WITH CLEAN WATER AND PLEASE TAKE IMMEDIATE ACTION AND CONSULT A DOCTOR.
- STOP IMMEDIATELY TO USE THE BATTERIES IN CASE

YOU NOTICE ANY ANOMALY AND AVOID THEIR CONTACT WITH METAL OBJECTS, SUCH AS KEYS OR COINS.

- NEVER DISASSEMBLE, CUT, COMPRESS, FOLD, PIERCE OR DAMAGE THE BATTERY. IN CASE YOU NOTICE SOME LIQUID LEAKING OUT FROM THE BATTERY AND AVOID IT TO CONTACT WITH YOUR SKIN OR EYES. IF IT HAPPENS, RINSE IMMEDIATELY THE AFFECTED AREA WITH ABUNDANT WATER OR CONTACT YOUR DOCTOR PROMPTLY. DO NOT ATTEMPT TO INSERT ANY FOREIGN OBJECT INTO THE BATTERY AND NEVER IMMERSE OR EXPOSE IT TO WATER OR ANY OTHER LIQUIDS.

- DO NOT USE BATTERIES AND THE BATTERY CHARGER FOR ANY OTHER USE EXCEPT ITS INTENDED ONE. AN IMPROPER USE OR THE USE OF NON-APPROVED OR NON-COMPATIBLE BATTERIES OR CHARGER MAY LEAD TO DANGER OF FIRE, DEFLAGRATION OR OTHER DANGERS FOR WHICH THE MANUFACTURER IS NOT RESPONSIBLE AND CANNOT PROVIDE ANY GUARANTEE.

- NEVER USE A DAMAGED BATTERY OR A DAMEGED CHARGER. USE THE CHARGER ONLY INDOORS.

- NEVER THROW BATTERIES INTO HOUSEHOLD WASTE.

- ALWAYS COMPLY WITH LOCAL REGULATIONS. RECYCLE WHENEVER POSSIBILE. DANGER OF DAMAGES AND / OR LEAKAGES!

- NON-COMPLIANCE WITH THE ABOVE-MENTIONED WARNINGS COULD CAUSE ELECTRIC SHOCK RISKS, DAMAGES / LEAKAGE OF LIQUIDS, OVERHEA-

TING, EXPLOSION OR FIRE, WHICH MAY LEAD TO SERIOUS PERSONAL INJURIES AND / OR DAMAGE TO OBJECTS.


DISPOSAL OF BATTERIES / ACCUMULATORS AT THE END OF LIFE

- THE BATTERY/ACCUMULATOR AND THE EVENTUAL ELECTRIC AND ELECTRONIC DEVICE MUST NOT BE DISPOSED AS URBAN WASTE BUT IN ACCORDANCE WITH 2002/96/CE, 2006/66/CE DIRECTIVES AND NATIONAL LAWS RELATING TO THE DISPOSAL OF THE PRODUCT AT THE END OF ITS LIFE. IT HAS TO RESPECT A SEPARATE COLLECTION IN ORDER TO AVOID ENVIRONMENTAL DAMAGES.
- IN ORDER TO REMOVE BATTERIES, PLEASE REFER TO THE SPECIFIC GUIDELINES.

NOTE:


FAILURE TO RESPECT WHAT STATED ABOVE IS SANCTIONED BY THE LAW.

1. PARTS DESCRIPTION

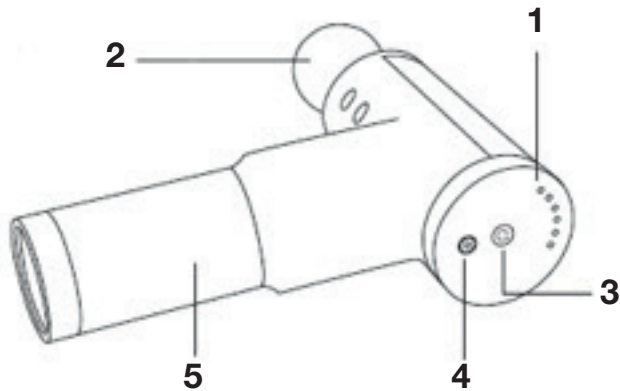1. LED AND FREQUENCY INDICATOR LIGHTS
2. REPLACEABLE MASSAGE HEAD
3. ON / OFF SWITCH
4. CHARGING INTERFACE
5. NON-SLIP GRIP

2. PRODUCT CHARACTERISTICS

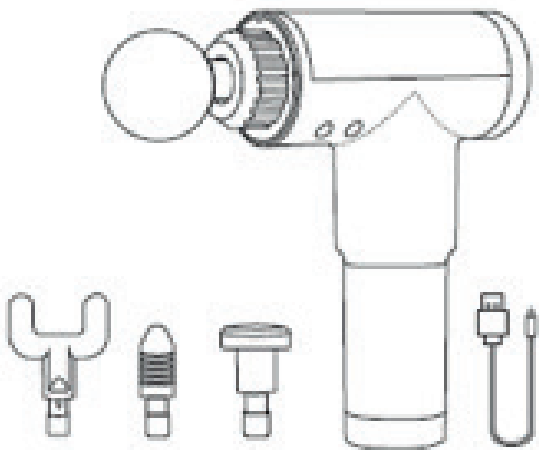6 ADJUSTABLE FREQUENCIES

4 REPLACEABLE HEADS FOR DIFFERENT KINDS OF MASSAGE

3. HOW TO USE THE MASSAGE GUN

Charging and Use

1. Place the device on a flat, stable and non-slip surface, then insert the connector of the recharging cable in the appropriate socket on the back of the device, then connect the plug of the power supply to the electrical outlet. Charge for about 2 hours before use.

2. To use the massage gun, disconnect the battery charger, insert the relevant massage head fully into the appropriate slot, without forcing, and press the button in the center of the rear to switch on the device. Successive presses of the switch allow you to select the percussion frequency. The lights on the circumference will illuminate to indicate the selected operating speed. To

turn the device off, press the center switch.

3. In order to prolong the life of the battery, it is recommended that you do not fully deplete the battery charge.


4. To remove the massage head, turn off the massage gun, and gently pull the attachment lengthwise. Use an antibacterial product or alcohol to clean the massage heads between uses.

WARNING!


- Do not use the device for purposes other than those intended!
- Do not use the product for more than 1 minute on the same area and in any case for intervals longer than 10 minutes in order to avoid overheating. Allow proper cooling. Wait at least 5 minutes before resuming use and starting a new maximum interval;

3. METHOD OF REPLACEMENT OF MASSAGE HEAD

1. Pull out the original massage head

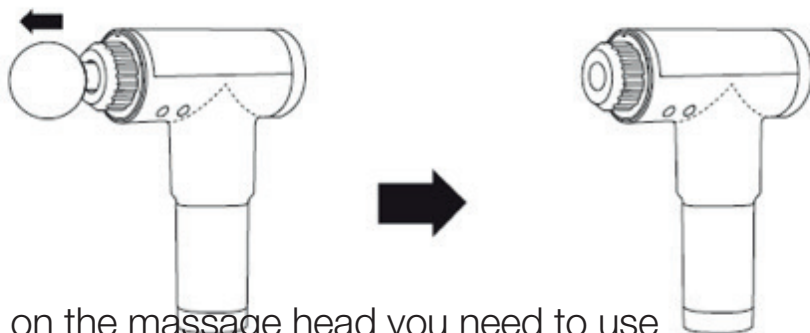2. Put on the massage head you need to use

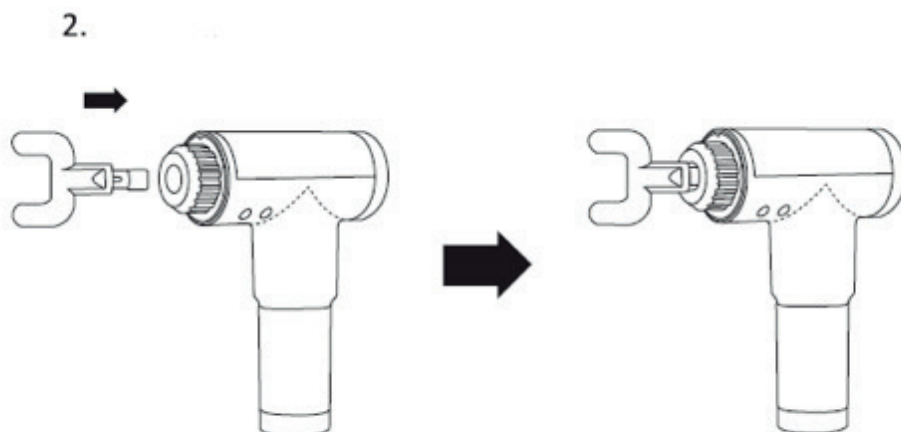4. WARNINGS

- For adult use only. Keep out of reach of children;
- Stop using it immediately in case of pain or discomfort;
- Avoid contact with the head and/or other bony areas;
- Do not use it for more than 1 minute on a single specific area of the body. Continue to move the gun across muscles and avoid bony areas;
- Avoid getting fingers, hands, and/or hair near the massage gun head;
- Do not drop the massage gun;

WARNING:

DO NOT USE THE MASSAGE GUN WITHOUT A DOCTOR'S APPROVAL UNDER THE FOLLOWING CONDITIONS:

- Recent surgery of any kind
 - Skin diseases and / or injuries caused by vibrating massage guns
 - Disc hernia
 - Pregnancy
 - Diabetes / neuropathies or retinal damage
 - Pacemakers
 - Wounds or injuries
-
- Do not immerse in water, sand or other loose debris that could enter the motor housing vent.
 - Do not attempt to disassemble the product or tamper with the product to alter its performance.

5. ADVERSE REACTION

Skin irritation and minor bruising may occur if the massage gun is used in a localized spot for more than 60 seconds or on bony areas. Adverse reactions may occur on persons recovering from surgery or persons with skin disease, neurological disorders, wounds, bone fractures and/or bruises. As a precautionary measure to avoid negative reactions, you should consult your physician before using the massage gun.

6. CLEANING

- To clean the device, use a slightly damp cloth, then dry thoroughly.
- Clean the heads using antibacterial products to prevent any bacterial contamination.

7. TROUBLESHOOTING

- If the massage gun does not turn on, check the charging connections, and proceed with a full charge.
- If the battery does not charge, check if the power cable is correctly connected to the charger or to the power outlet.

8. BATTERY AND ADAPTER SPECIFICATIONS

Battery specifications:

2 x 18650 cells,

Voltage: 7.4 V

Wattage: 10 W

Battery capacity: 1500 mAh

Charging time: about 2 h

Battery life: about 40 minutes

Adapter specifications:

Input 100-240 V AC | 50/60 Hz | 0.4 A

Output: 8.4 V DC | 1.0 A

Note on battery life:

If the device is not used for a long period of time, it is recommended to charge the battery to 50% every 3 months.